




Not Food and Drug Administration (FDA) approved for this use For Forensic Use Only

QUICK REFERENCE INSTRUCTIONS

Step 1

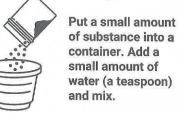

The amount should cover the president's hair on a penny.




Step 2

Dip the end of the strip for 15 seconds.

Lay the strip flat.

Step 3



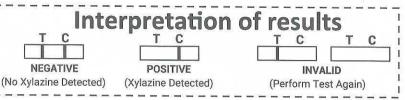

WARNING: The results of the SAFe-X test are not a guarantee that the tested substance can be safely consumed. Even if the results of the SAFe-X test are negative for xylazine, anyone consuming the tested substance could overdose, sustain severe injury, or die. Failure to follow the above instructions could result in false positive or false negative test results. The above instructions are only for the SAFe-X test, and you should not use them for any tests other than the SAFe-X test.

Please see the reverse side for other warnings and information regarding the SAFe-X test.

Please see other side for Quick Reference Instructions





Not Food and Drug Administration (FDA) approved for this use For Forensic Use Only

Read, understand, and follow the instructions on the reverse side your fa

- Read, understand, and follow the instructions on the reverse side, your failure to do so could result in false positive or false negative test results, severe injury or death.
- Even if the results of the SAFe-X test are negative for xylazine, anyone consuming the tested substance could overdose, sustain severe injury, or die.
- Do not consume illicit drugs, especially if you suspect that they contain, or the SAFe-X test result is positive for, xylazine.
- The SAFe-X test does not detect the presence of substances other than xylazine, and it does not detect all forms of xylazine.
- The SAFe-X test does not detect the percentage of xylazine present in the tested substance.
- Exposure to even small amounts of xylazine, or handling of xylazine, can cause severe injury or death.
- When testing a substance, use protective equipment, such as gloves and eye protection.
- Do not use the SAFe-X test after expiration or if the packaging is damaged or shows other signs of tampering.
- Dispose of the SAFe-X test in accordance with applicable law.
- Anyone using the SAFe-X test accepts all responsibility for any injury, death, criminal penalty, or civil liability that occurs from
 consumption of a substance, whether an illicit drug or not, and regardless of whether the substance was tested, or was not
 tested, for xylazine.
- Call 911 or seek immediate medical attention if you think you or someone else has overdosed or is at risk of overdose. Common signs of an overdose are listed in the "How do I recognize an overdose?" section below.
- If you or someone you know requires counseling or treatment for drug use, call the SAMSHA National Helpline at 1(800)-622-HELP(4357). This helpline is confidential and available 24/7, 365 days a year.

How do I recognize an overdose?

Signs of an overdose include:

- Small, constricted "pinpoint pupils"
- · Falling asleep or loss of consciousness. Slow, shallow breathing
- Choking or gurgling sounds
- Limp body
- Pale, blue, or cold skin
- Because this list is not exhaustive, call 911 or seek immediate medication attention if you think you or someone else has overdosed or is at risk of overdose.